




ALL THE ESSENTIALS FOR EYE EXAM

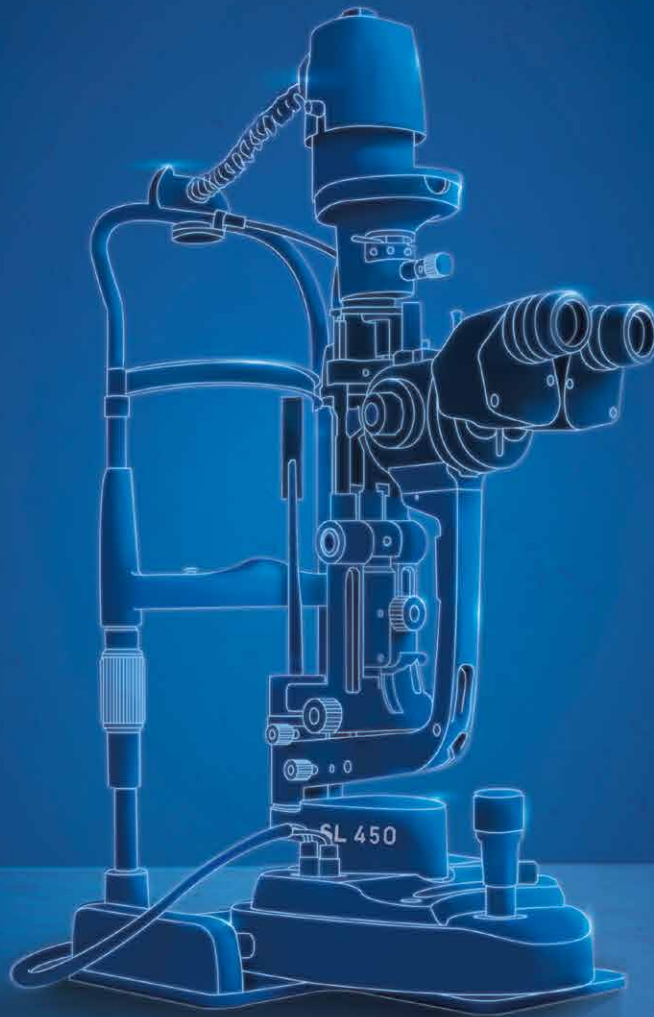

SL 350 / SL 450

SLIT LAMPS





HD OPTICAL SYSTEM

- High definition optical resolution providing a wider and better vision for more details of the pathologies.
- Magnification levels easily accessible by drum rotation.
- Instant access to all functions thanks to its compact, ergonomic design.



BRIGHT LED ILLUMINATION

- Tower illumination systems, featuring up to 20° tilting tower illumination.
- Low color temperature and high illumination intensity for better experience.
- Longer life span and less energy.



BUILT-IN YELLOW FILTER

- Integrated yellow filter along with cobalt-blue filter increases the contrast of Fluorescein Staining image for your daily use.



EASY INSTALLATION

- Compact control panel with built-in power supply makes it much easier to install while keeping a soft practice.

OPTIONAL ACCESSORIES

- Measuring eyepiece, observation tube, tonometer adaptor.

SPECIFICATIONS

MAGNIFICATION	
Type	Galilean
Cylinder	12.5x
Selection	Drum rotation
Field of view	
SL350 L	10x (22mm) - 16x (14mm) - 25x (8.9mm) 6.3X (36.2mm) - 10x (22mm)
SL450 L	16x (14mm) - 25x (8.9mm) 40x (5.7mm)
BASE	
Longitudinal movement	110 mm
Lateral movement	110 mm
Vertical movement	30 mm
Chinrest (vertical)	80 mm

SLIT	
Slit width	Continuous from 0 to 14mm
Slit length	Continuous from 1 to 14mm
Aperture	0.2 to 14mm
Slit angle	0 to 180°
Tilt	5, 10, 15, 20°

ILLUMINATION LAMP	
SL350L/SL450L	3V,3W LED lamp

FILTERS - Heat-absorbing, ND, Red-free, Cobalt blue, Built-in yellow filter	
Power-supply	220/110V (±10%)
Dimension and weight	740mm x 450 mm x 530 mm (L/W/H) 22kg

As improvements are made, these specifications and pictures are not contractually binding and may be changed without prior notice.



ESSILOR INSTRUMENTS
81 boulevard Jean-Baptiste Oudry
94000 Créteil
France
Tel. : +33 (0)1 49 80 62 80
www.essilor-instruments.com

